

i05000


DENTAL LIGHT



 Belmont®


"Clean Station" to surgeries around world.

CLESTA iO 5000

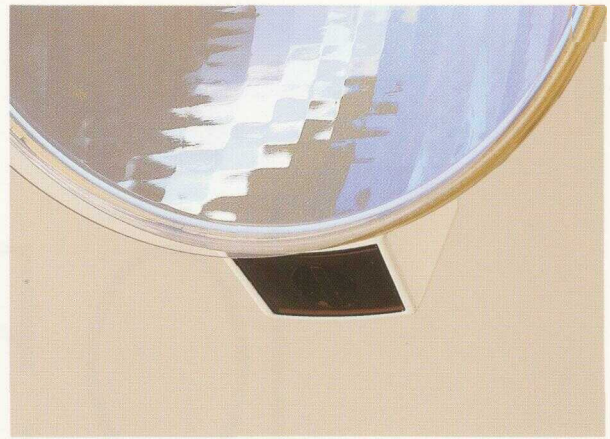

A movement of the hand, beneath the switch, turns the light on and off, offering protection to the dentist and assistant from the risk of cross-contamination.

Approximating the ideal white of natural daylight. Color shade matching, with a light color in the 5000° K range, will offer a more accurately perceived color of an object. Unique sensor switch allows touchless on and off of the light.

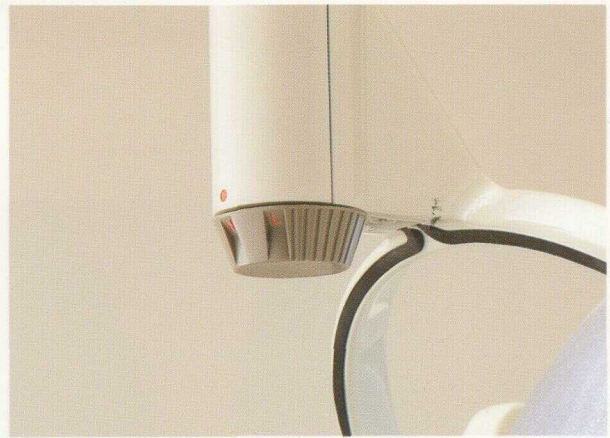

Three-position controller located in the rear of the light head enables the selection of three intensities, low, medium or high which eliminate any variations that might cause eye fatigue.



Smooth, balanced movement. Lightweight, long arms provide for easy positioning and greater reach. iO 5000 can be mounted on the unit or ceiling. Track light system is also available.

Specifications:

Rating Voltage: AC. 100/110/115V and 220/230/240V

Power Output: 65W

Output Voltage: 23.7V 21.9V 19.6V

Bulb Type: Tungsten Halogen Type (JA24V60W)

Light Output: High 28,000 lux

Medium 22,000 lux

Low 15,000 lux

Color Temperature: 5,000 Kelvin

Light Pattern: 70mm x 190mm at 750mm

The manufacturers reserve the right to alter and improve the specification without notice.

 **Belmont**[®]

Printed in Japan 9912MA